


# BIORING®

ADJUSTABLE GASTRIC BANDING  
FOR MORBID OBESITY TREATMENT

## Product codes and sizes :

	Product code	Size (cm)	Description	Packaging
	RING 230 000	Ø 2,7 / 2,9	+ Syringe + Catheter rincer	1 / box
	RING S20 000	Ø 2,1 / 2,4	+ Gripping tab	1 / box
	RING PORT P1	Small	Polysulfone port + Connection system	1 / box
	RING PORT P3	Large	+ Huber needle	1 / box
	RING PORT T1	Small	+ <b>Titanium</b> port + Connection system	1 / box
	RING PORT T3	Large	+ Huber needle	1 / box
	Product code		Description	Packaging
	RING CATH 00		Stomach tube	1 / box
	RING KIT REM		Replacing kit	1 / box

Non contractual document. Specification and photos can be modified without informations. docBIORING.GB - 03/04/07  
Cousin Biotech S.A.S. with a capital of 320 000 € - 398 460 261 RCS Lille - N°TAT FR 34 398 460 261

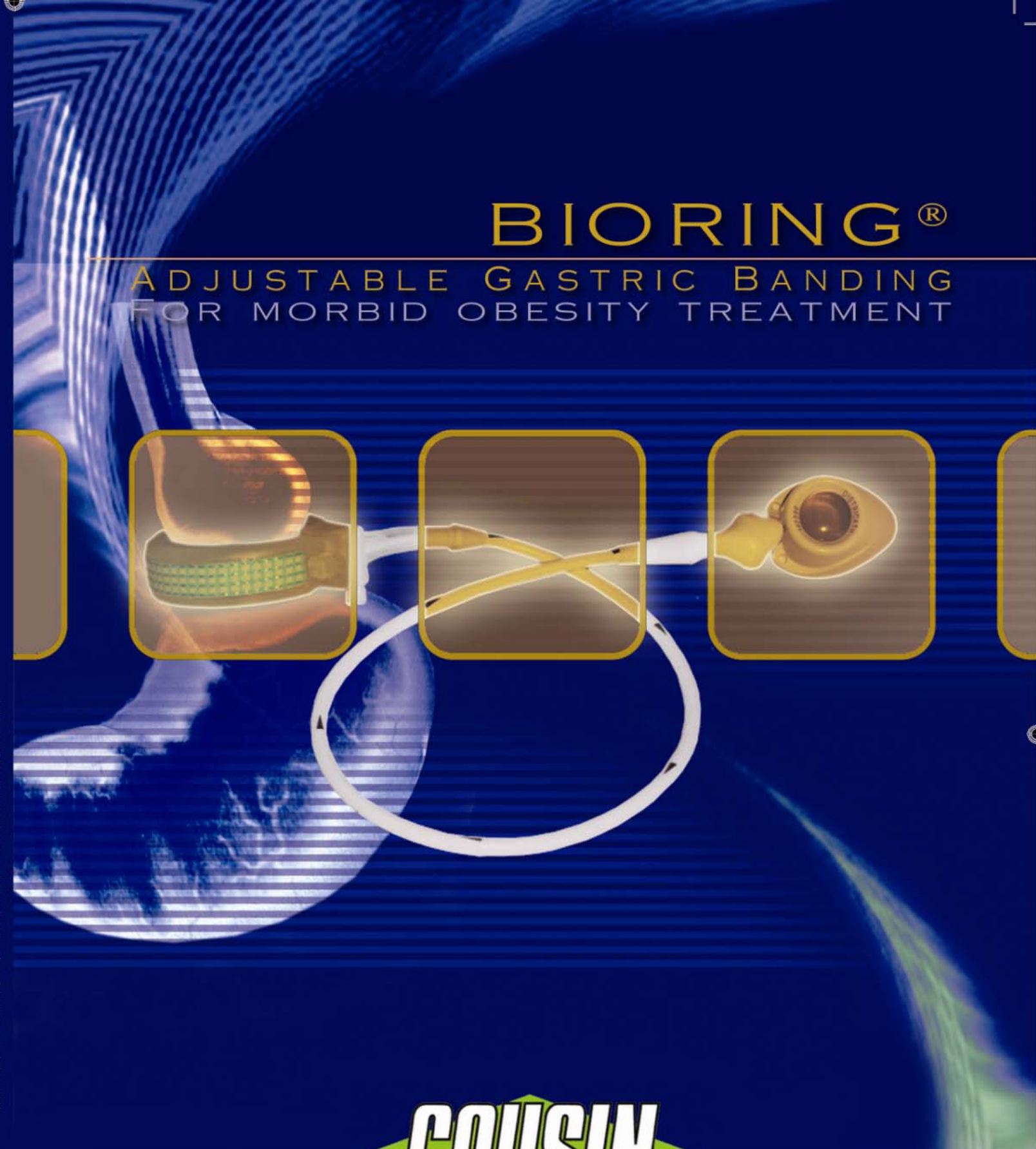

**COUSIN**  
B I O T E C H

Distributed by :

**COUSIN**  
B I O T E C H

COUSIN BIOTECH  
8, rue de l'Abbé Bonpain  
59117 Wervicq-Sud  
FRANCE  
Tel. : +33 (0) 3 20 14 40 00  
Fax : +33 (0) 3 20 14 41 22  
mail : visceral@cousin-biotech.com  
www.cousin-biotech.com




# BIORING®

## ADJUSTABLE GASTRIC BANDING FOR MORBID OBESITY TREATMENT



Preshaped body.

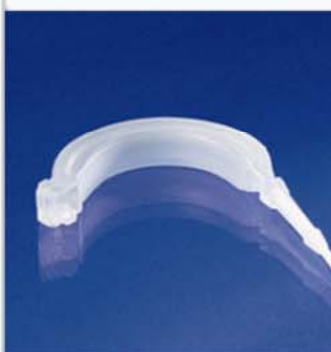

Low pressure expandable system.



X-rays detectable ring and catheter.



Marked catheter.



4 sizes of implantable ports.

### ➤ FLEXIBLE GASTRIC BANDING

Dimethyl siloxane (Silicone).  
Low pressure filling system.

### ➤ EXPANDABLE SYSTEM

Avoids reversal  
and slippage of the ring.

### ➤ PRE-SHAPED BODY

Easy surrounding.

### ➤ LOW PRESSURE BALLOON

Covering perfectly the gastric wall  
after placement.  
Avoids risks of nipping.

### ➤ LENGTH OF THE BODY

2 sizes: 12cm or 10,5 cm, both are adapted to  
every kind of surgical technique such as  
pars flaccida or peri gastric procedure.

### ➤ WIDE VESTIBULE

Once inflated the contact surface increases  
from 10,6 cm<sup>2</sup> to 13 cm<sup>2</sup> for the large ring,  
and from 9,3 cm<sup>2</sup> to 11,4 cm<sup>2</sup> for the small one,  
allowing an optimal pressure repartition on the  
gastric wall and decreasing the erosion risk.

### ➤ TROCAR

12 mm.

### ➤ CLOSURE SYSTEM

Easy and efficient.

### ➤ POSSIBLE RE-OPENING

For a perfect positioning.

### ➤ NO FLASH

Better visibility.

### ➤ X-RAYS DETECTABLE RING AND CATHETER

Easy radiographic control.

### ➤ MARKED CATHETER

Makes easy its orientation.

### ➤ LOCKING SYSTEM

Secured connexion between  
port and catheter.

### ➤ IMPLANTABLE PORTS

4 sizes