


Siesta *i* Whispa

The Integrated Solution



*i*ntegrated Breathing System
*i*nnovative user friendly design
*i*nspiring Anaesthesia Workstation

www.dameca.com

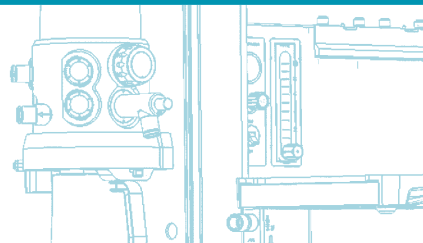

Siesta *i* Whispa



- New Whispa ventilator
- Integrated **B**reathing **S**ystem (IBS)
- Suitable for low flow/minimal flow
- Minoxygard (min. 25% O₂)
- Low maintenance costs
- Low lifetime costs
- Easy attachment of patient monitor
- Ceiling pendants/wall versions available

- 8" colour TFT screen or monochrome 7" LCD screen
- Display with graph and loop
- 5 ventilation modes: VCV, PCV, SIMV, manual & spontaneous.
- Microprocessor controlled ventilator
- Spirometry
- Built-in oxygen analyzer
- Same bag-in-bottle for adults and infants
- Automatic self test



Features



The **Minoxiguard** function ensures a minimum 25% oxygen in O₂/N₂O mix, and is further equipped with a slave regulator, cutting off N₂O and giving audible alarm in case of O₂ supply failure. Minoxiguard is also available with twin/dual flow tubes for O₂ and N₂O.



The new **IBS (Integrated Breathing System)** is a fusion of the bag-in-bottle, absorber, and patient circuit. These 3 parts are now integrated into one system, thus drastically reducing tubing connections between machine and patient. Built-in manual ventilation valve (5-50 hPa), just above the manual ventilation bag connection.



Siesta i Whispa is equipped with various safety/emergency features, such as fresh gas flow switch between circle system (IBS) and auxiliary gas outlet. Additional features: O₂ flowmeter with separate outlet with tube nipple and integrated patient suction unit.



Back-bar for 2 selectatec vaporizers (German or British standard) with interlock safety mechanism.

Dameca



The Siesta *i* Whispa is an anaesthesia workstation, developed and manufactured by Dameca a/s, Denmark, a leading independent Scandinavian manufacturer of anaesthesia workstations and equipment. Established in 1947, Dameca a/s has concentrated its development and manufacturing activities solely on anaesthesia machines, ventilators and accessories for the Danish and world market. Over the years, Dameca a/s has built up extensive knowledge regarding precision manufacturing - here illustrated out latest production center with robot. We are very proud that many anaesthetists world-wide pass on to us their valuable experience and requirements to ensure that Dameca a/s can always offer the right anaesthesia solutions.



Please contact us or our local distributor for further information regarding Siesta *i* Whispa or Dameca's other anaesthesia machines.

We are constantly developing our products and consequently reserve the right to make alterations in design and equipment without prior notice. 11/2004